
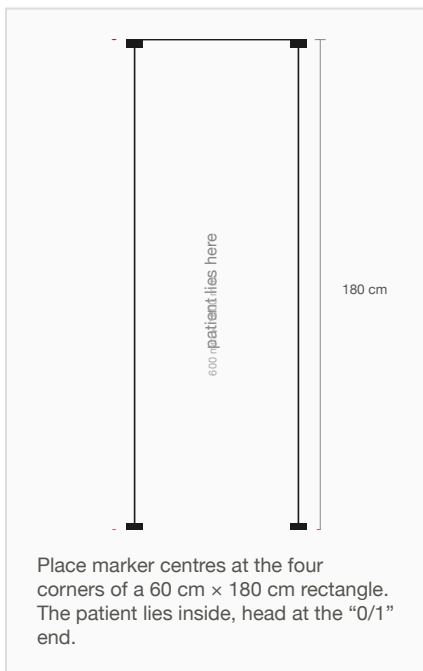
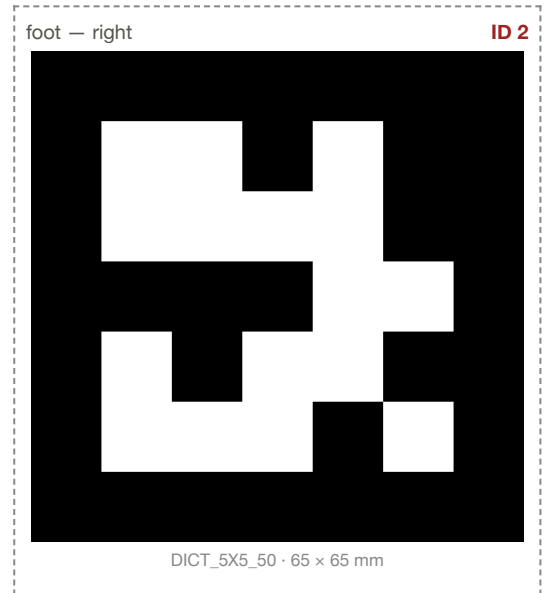
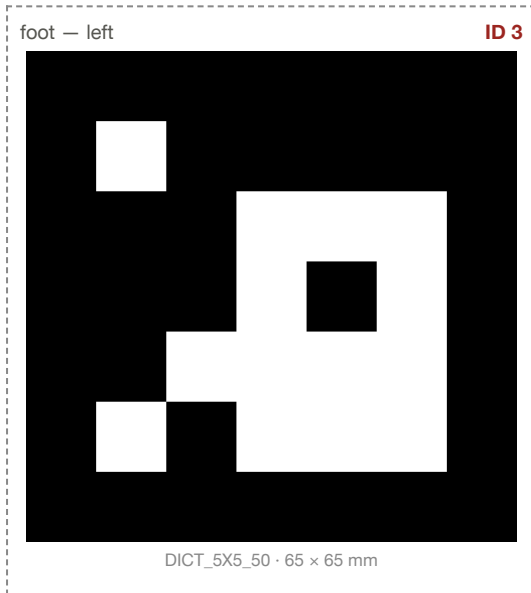
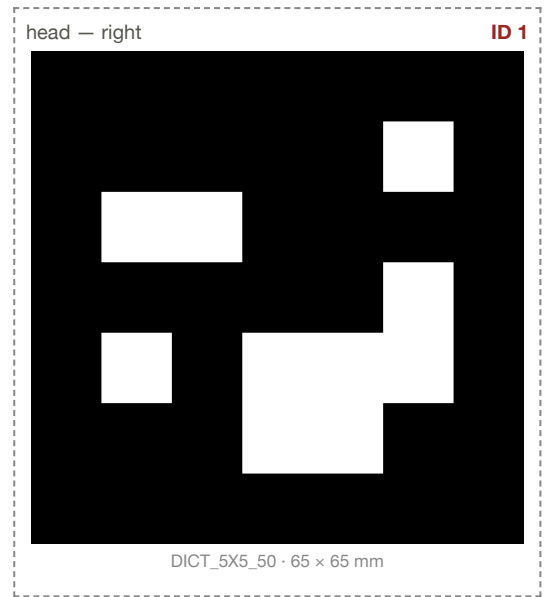


# Bed-Mat Marker Sheet

Print at 100% scale (no “fit to page”), cut along dashed lines, place on bed/floor as shown.

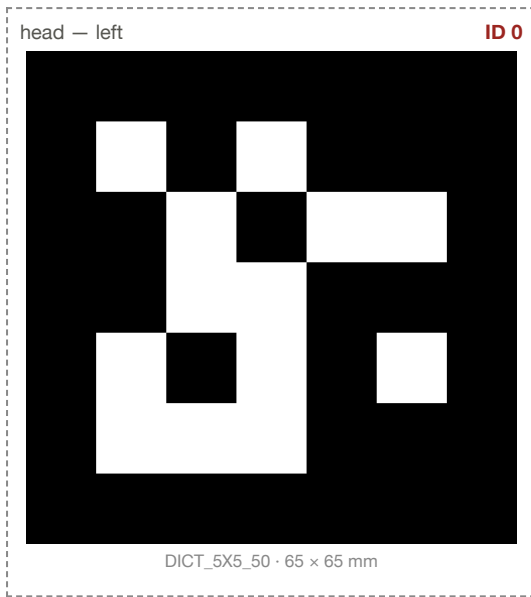

## Setup steps

1. Print this sheet at **100 % scale**. Verify with the ruler check below — the bar should measure exactly 5 cm.
2. Cut along the four dashed boxes.
3. Lay the four markers face-up on the bed mat or floor in the rectangle shown to the left, IDs in the labelled positions.
4. Patient lies supine inside the rectangle, head at the 0/1 end, arms by their sides.
5. Take one overhead photo so all four markers are clearly visible — phone on a chair / tripod aimed straight down works fine.
6. Upload to the anthroheight UI; with all four markers detected the system computes the bed-plane scale automatically.

**Scale check — this bar should measure exactly 5 cm**

